

LECKEY

Product Range Guide



Say hello to the
Leckey family of products.



I am delighted to introduce our new Leckey product guide and celebrate our Leckey family circle.

At Leckey we believe in **'Supporting Children'**.

To everyone in Leckey, Supporting Children means working together with our family of experts – designers, therapists, parents and the broader clinical community to achieve our goals.

Our **Mission** is; To improve the quality of life of children with special needs.

Our **Vision** is; A world where each child is included and they can access their full potential.

Our **Courageous Goal** is; To prevent hip migration in children with CP by 2025.

In Leckey we see hip migration as the single biggest life event that our children may experience and has the greatest impact on their quality of life and their carers from that moment on.

Everything we do in Leckey, in partnership with our family of experts, is aimed at achieving these goals.

Our approach is simple – lead through innovation and clinical focus.

Through early intervention, childhood right up to adulthood, Leckey's experienced team of designers and clinicians work together with our family of experts to develop exciting innovative products that make a real difference to the quality of life for the children and their families.

We aim to develop clinically and technically superior solutions, products and programs, that improve the child's Mobility, Ability and Participation.

All of our designs undergo extensive clinical and usability trials with physiotherapists, occupational therapists, users and their families.

They are supported by in-depth clinical research studies where we work with leading universities and clinicians internationally to ensure that all Leckey products deliver on the clinical needs of the healthcare professionals and the social needs of the user.

We are committed to our Mission, Vision and Courageous Goals. We want to work with you in partnership to build our family - a community of experts with shared goals and values.

Join our family and together we will make a difference.

- The Leckey Team

Contents



Seating and Mobility

05

- BeMe
- Early Activity System
- Pal Seat
- Squiggles Seat
- Therapy Bench



Standing and Walking

17

- Squiggles+
- Squiggles TT
- Mygo Stander
- Horizon
- Pronestander
- Freestander
- Totstander
- MyWay+
- MyWay Pedal



Other Products

33

- Advanced Bath Chair
- Shower Trolley
- Sleepform

A thick white arc that starts from the left edge of the page, curves upwards and then downwards, ending at the right edge. It is positioned above the main title.

Seating and Mobility

LECKEY BeMe

The infinitely configurable BeMe, is a revolution in supportive seating that puts a child's development and experience at the heart of its design. Empowering a child to simply BeMe.





**BeMe delivers
posture, function
and comfort without
compromise**

As Unique As A Fingerprint

Building on decades of Leckey clinical and design expertise and responding to detailed feedback from therapists and carers, BeMe is a highly personalised seating system, available in three sizes, and spec'd to perfection. Unique modularity provides easy configurability to fit the exact needs of the child as they change and grow. Promoting inclusivity in all day-to-day activities through good posture, BeMe encourages confidence and independence - empowering a child's beaming smile to shine through.

A Revolution in Seating

BeMe is a modular seating system comprising three sizes with three seat options: contoured, basic and split and three back supports: basic, moderate and complex to suit every need.

**Click and Switch
with Ease**

**Total ease of use and
compatibility are at
the heart of BeMe's
responsive design**



It's surprisingly simple to find the perfect fit with BeMe. It is cleverly designed so that you can configure the right support quickly and easily.

- Configuring BeMe is a snap as parts can be easily and intuitively interchanged
- Clinical adjustments, including seat depth and back height, can be made while the child is seated
- Armrest, femoral guides and pommels can be quickly released making transfers easier than ever

Revolutionary seating with world's first truly growable contoured cushion

BeMe's child-centred ethos inspired Leckey to create the first genuinely growable contoured cushion. This innovative growable cushion,

with specially developed high-stretch foam and fabric, is contoured to maximise pelvic stability and comfort as the child changes, grows and develops.

This breakthrough in support technology envelops the child for immersive comfortable positioning, improving pressure distribution and enhancing function.

By reducing the need for hardware like pommels and hip guides, BeMe allows the child to take centre stage without compromising support.



Product Features

- Manoeuvrable Hi-Low Mobility base that enhances participation and inclusion
- Modular seating to fit the exact needs of the child as they change, develop and grow
- Genuinely growable contoured cushion – providing support, skin protection and comfort even while child is growing and changing.
- Diverse range of accessories and options to provide the required postural support
- Optional air-controlled dynamic backrest to better manage a child's dynamic movement
- Now available on the Zippie Iris and the Zippie Q300M Mini



Zippie Iris



Zippie Q300M Mini

Available in 4 fabric colours

Grey and Blue options are stocked locally



Red and Purple options are not stocked locally and extended lead times may apply

Technical Details					
Seat	Small	Medium	Large	Medium Long	Large Long
Seat Depth	170*-325mm	275-390mm	350-475mm	325-440mm	400-525mm
Seat Width	275mm	325mm	415mm	325mm	415mm
Max User Weight**	35kg	50kg	70kg	50kg	70kg
Backrest	Small	Medium	Large		
Back Height	360-470mm	390-520mm	500-650mm		
Back Height with Contoured Seat	370-480mm	400-530mm	475-625mm		
Lateral Width	130-230mm	200-300mm	225-350mm		
Hip Guide Width	165-230mm	200-260mm	250-350mm		
Backrest Angle	+10° to -25°	+10° to -25°	+10° to -25°		
Max User Weight **	35kg	50kg	70kg		
Footrest	Small	Medium	Large		
Lower Leg length with Basic/Split Seat	185-285mm	280-450mm	290-510mm		
Lower Leg length with Contoured Seat	185-285mm	290-460mm	300-520mm		
Maximum Dorsiflexion/Plantarflexion	+/- 10°	+/- 10°	+/- 10°		
Chassis	Small	Medium	Large		
Maximum Weight Limit**	35kg	50kg	70kg		

*170mm-220mm seat depth required spacer options

**System max user weight should be lesser of seat, chassis and backrest

LECKEY Early Activity System

Colourful, comfortable, supportive activity system for children aged 0-48 months.



Product Features

- Offers a gentle, playful introduction to 24-hour postural management.
- Practical all-in-one kit, providing versatile support. The positioning possibilities are endless and everything can be packaged up in the duffle bag for easy transport.
- Supplied with the popular Activity Development Programme which supports parent-led therapy at home.
- Soft-touch, anti-bacterial, machine washable and tumble dry safe covers.



Key Positions

- Back Lying (Supine)
- Tummy Lying (Prone)
- Side Lying
- Floor Sitting (Long sitting)
- Hands and knees (Four point kneeling)

Technical Details		
One Size Only		
Age	Max 36 months	
Weight	Max User Weight 15kg / 33lbs	
Stature	Max 1000mm / 39.4 inches	
Floor Mat	Width	"760mm / 26.3 inches"
	Length	950mm / 37.4 inches"
Sitting Support	Height	"190mm / 7.4 inches"
	Depth	300mm / 11.8 inches"
Seat width	Min	"200mm / 7.9 inches"
	Max	250mm / 9.8 inches"
Support Roll small	Width	"60mm / 2.4 inches"
	Length	600mm / 23.6 inches"
Support Roll medium	Diameter	"80mm / 3.2 inches"
	Length	700mm / 27.5 inches"
Trunk Wedge	Width	"270mm / 10.6 inches"
	Length	208mm / 8.2 inches"
	Angle	7.5 degrees"
Positioning Straps	Width	"80mm / 3.14 inches"
	Length	660mm / 25.9 inches"



LECKEY Pal Seat

Colourful, comfortable classroom seating system providing mild support for children aged 1-12 years.

Product Features

- Two-point lap belt and side pads for lateral stability. An optional ramped cushion can be added if extra help is needed to keep the pelvis and femurs in alignment.
- Where a less upright angle is required, the backrest angle and armrests adjust to suit the child's needs, and along with the large tray or table, provide appropriate support.
- Height adjustable armrests, legs and footplate.
- Available in two bright colours in wipe clean fabric [please see website for details].
- A range of base and footrest options available, including castors, rockers and stabilisers.



Available in 2 fabric colours

Technical Details				
	Size 1	Size 2	Size 3	Size 4
Age	1-3	2-5	4-8	7-12
User Height	Min 700mm [27.6"] Max 950mm [37.4"]	Min 900mm [35.4"] Max 1050mm [41.3"]	Min 1000mm [39.4"] Max 1300mm [51.2"]	Min 1200mm [47"] Max 1600mm [63"]
Max. User Weight	30kg	35kg	40kg	50kg
Seat Depth	180mm-240mm [7"-9.4"]	220mm-290mm [8.6"-11.4"]	250mm-330mm [9.8"-12.9"]	300mm-400mm [11.8"-15.7"]
Seat Width	150mm-250mm [5.9"-9.8"]	170mm-270mm [6.6"-10.6"]	230mm-330mm [9"-12.9"]	280mm-380mm [11"-14.9"]
Seat Height	160mm-240mm [6.2"-9.4"]	185mm-285mm [7.2"-11.2"]	250mm-350mm [9.8"-13.7"]	300mm-400mm [11.8"-15.7"]

LECKEY Squiggles Seat

Suitable for use as a home chair, school chair or with a Zippie base, Squiggles meets all the postural and comfort requirements of a seating system for young children with special needs, while maintaining a fun, colourful design.

Product Features

- Designed specifically for the developmental needs of younger children- the one-piece seat base and adjustable sacral support can be used to accommodate flexion or encourage extension.
- The shape of the seat and back, the laterals with their fabric tip which wrap securely around the child and the way in which the 4-point pelvic harness integrates with the hip pads, provides a snug, comfortable fit and streamlined look.
- Being low to the ground, the Squiggles Seat promotes peer to peer interaction, making it ideal for home and nursery/school environments.
- Adjustable footplate can be tilted to the preferred angle, ensuring the feet are well supported and providing a secure base for maximum upper body function. The footplate flips away to allow for easy transfers.
- Tilt-in-space feature provides a simple way to achieve a change of position without disturbing overall posture.



Available in 3 fabric colours

Technical Details	
One Size Only	
Age	1-5
User Standing Height	750mm-1100mm [29.5"-43.3"]
Max. User Weight	22kg [48.4lbs]
Seat Depth	190mm-290mm [7.5"-11.4"]
Seat Width	160mm-260mm [6.3"-10.2"]
Chest Width	160mm-220mm [6.3"-8.75"]
Backrest Height	320mm-380mm [12.6"-15"]
Lower Leg Length	122mm-270mm [4.8"-10.6"]
Tray Height	115mm-170mm [4.5"-6.7"]
Backrest Angle	Prone 10° - Recline 25°
Seat Height (Hi-Low Chassis)	270mm-650mm [10"-25"]

**Squiggles Seat now able to
interface with Zippie IRIS
Crash Tested Combination!**



LECKEY Therapy Bench

Developed to support children in a neutral posture, this therapy bench with innovative activity cradles facilitates child active therapy by freeing up the hands.



Product Features

- Innovative, corseted activity cradle stabilises the pelvis in a neutral position, acting as an extra pair of hands for the therapist whilst giving the child the confidence to move, reach, and grasp.
- Height adjustable trunk support helps the child to sit upright with the head and spine stacked over the pelvis. This is the best and easiest way to learn sitting posture.
- Supportive, height adjustable therapist pad acts as extra support for the therapist who can rest against the rear support and gently encourage shoulder retraction.
- Facilitates lateral and anterior/posterior tilt to promote child activity.
- Height and depth adjustable bench making it versatile for both larger and smaller users.
- Acts as an aid to therapists in the Segmental Assessment of Trunk control (SATCo), providing a greater understanding of where to apply postural support to optimise function.



Therapy Bench

Technical Details	
One Size Only	
Max User Weight	80kg [176.3lbs]
Age	1-14
Bench Footprint	620mm x 740mm [63.7" x 29.1"]
Seat Base Height (from ground)	Min 330mm [12.9"] Max 460mm [18.1"]
Seat Base Width	620mm [24.4"]
Seat Base Depth	360mm [14.1"]
Footplate Height (from ground)	Min 20mm [0.7"] Max 110mm [4.3"]
Seat Base to Footplate Height	Min 220mm [8.6"] Max 440mm [17.3"]

Activity Cradle

Technical Details		
Cradle Size	Hip Circumference (mm)	ASIS - ASIS (mm)
Size 1	Min 400mm [15.7"] Max 500mm [19.6"]	Under 120mm [4.7"]
Size 2	Min 440mm [17.3"] Max 575mm [22.6"]	Min 130mm [5.1"] Max 160mm [6.2"]
Size 3	Min 520mm [20.4"] Max 650mm [25.5"]	Min 160mm [6.2"] Max 180mm [7.0"]
Size 4	Min 600mm [23.6"] Max 740mm [29.1"]	Min 190mm [7.4"] Max 200mm [7.8"]
Size 5	Min 660mm [25.9"] Max 855mm [33.6"]	Over 200mm [7.8"]





LECKEY



Standing and Walking

LECKEY Squiggles+

The Squiggles+ retains the best of the Squiggles stander but also offers heaps of added benefits so little ones can continue to enjoy standing and be part of the action.



Product Features

- Anatomically-correct abduction to stabilise the joint and provide optimum loading for early hip joint development. Squiggles+ offers 60° of hip abduction bilaterally in prone, upright or supine. It promotes proper alignment in the hips, knees and feet as pivot points are in line with the hips.
- Increased comfort and support with wrap-around knee supports. The comfortable medial knee supports distribute pressure and are especially suited to abducted standing. The padded wraps maintain the knee in neutral alignment.
- Providing greater confidence with transfers- the larger hip pads improve the carer's and child's confidence and feeling of security when transferring. The improved pivot chassis facilitates full horizontal loading.
- Activity tray with improved height and adjustability.
- Height-adjustable, foldable push handle.



Technical Details

	One Size Only
Age	1-5
Max User Weight	22kg/48.4lbs
User Height	Min 730mm / 28.5 inches Max 1110mm / 43.5 inches
Chest Support	Angle adjustment +/- 30° Depth adjustment 50mm / 2 inches
Chest Support Height (footplate to armpit)	Min 560mm / 22 inches Max 940mm / 37 inches
Chest Width (distance between laterals)	Min 160mm / 6.3 inches Max 230mm / 9 inches
Hip Support (height from footplate to mid buttocks)	Min 300mm / 11.8 inches Max 530mm / 20.8 inches
Hip Width (distance between hip guides)	Min 160mm / 6.3 inches Max 230mm / 9 inches
Distance between Midline of Knees	Min 140mm / 5.5 inches Max 215mm / 8.4 inches
Footplate Angle	Plantarflexion 10° Dorsiflexion 10°

Available in 3 fabric colours and a vinyl version



LECKEY Squiggles TT

Evolutionary early intervention standing frame that facilitates Targeted Training for children aged 1-5 years.



Product Features

- Maintaining the great features of the Squiggles+, the Squiggles TT goes a step further and facilitates Targeted Training. This evidence-based therapy is carried out with the child in an upright position, training one segment of their trunk at a time, allowing them to gain sequential control from the top down. Improved trunk control will help increase functional skills e.g. self-care and play. Better head control can help stop a child dribbling, improve their eye control as well as their ability to communicate.
- Squiggles TT is the only product specifically designed for Targeted Training therapy and which caters for all children. Horizontal loading is possible for complex children with poor head control and Squiggles TT is the only product to target all 7 trunk segments.
- Offers a fun and easy way for parents to carry out therapy at home. The snug supportive wraps are easy to use, the ratchet system with repeatable locking eliminates the fear of changing settings and the play ideas booklet provides support for home use.
- Dynamic base allows children to experience rocking movement whilst gaining upright postural control.
- Padded tray facilitates Targeted Training at the upper thoracic level.
- Conversion kit available for existing Squiggles+ users.



Available in 3 fabric colours and a vinyl version

Technical Details	
One Size Only	
Age (approx.)	1-5
Max User Weight	22kg/48.4lbs
User Height	Min 730mm / 28.5 inches Max 1100mm / 43.3 inches
Trunk Support Height (Distance from footplate to armpit)	Min 160mm / 6.25 inches Max 820mm / 32 inches
Trunk / Pelvic Width	Min 135mm / 6.25 inches Max 250mm / 9 inches
Pelvic Support Height (Distance from footplate to mid-buttock)	Min 300mm / 11.8 inches Max 530mm / 20.8 inches
Distance between midline of knees	Min 140mm / 5.5 inches Max 215mm / 8.4 inches
Hip Width (distance between hip guides)	Min 160mm / 6.3 inches Max 230mm / 9 inches
Distance between Midline of Knees	Min 140mm / 5.5 inches Max 215mm / 8.4 inches
Footplate	Plantarflexion -10° Dorsiflexion -10°
Pivot chassis angle range	90-180°
Tray Height	Min 480mm / 18.5 inches Max 825mm / 32.5 inches

LECKEY Mygo Stander

Strong, practical and adjustable standing system for kids aged 4-14 years.



Product Features

- 3-in-1 prone, supine and upright stander.
- Designed to improve upper limb function. The chest pad design replicates the shape of the rib cage which allows the shoulder girdle to freely protract and retract in prone, whilst in supine, the shoulder support promotes shoulder protraction. Furthermore, Mygo offers a lower tray height than most standers.
- Accommodates hip and knee flexion contractures of up to 25 degrees. The movement of the footplate and knee supports are linked which allows the user to stand with hips and knees flexed. This results in a more upright posture combined with the ability to support users whose contractures previously made them unsuitable for standing therapy.
- Use of the pelvic derotation belt puts the user's pelvis in a central and neutral position with minimum effort.



Available in 3 fabric colours

Technical Details		
	Size 1	Size 2
Age	4-10	8-14
User Height	1080mm-1450mm [42.5"-57"]	1350mm-1700mm [53.1"-66.9"]
Max. User Weight	50kg (110lbs)	60kg (132lbs)
Chest Support Height (footplate to armpit)	760mm-1080mm [29.9"-42.5"]	980mm-1300mm [38.5"-51.1"]
Chest Width (distance between laterals)	160mm-280mm [6.2"-11"]	200mm-320mm [7.8"-12.5"]
Hip Support Height (footplate to mid-buttocks)	460mm-730mm [18.1"-28.7"]	620mm-875mm [24.4"-34.4"]
Hip Width (distance between hip guides)	220mm-330mm [8.6"-12.9"]	280mm-390mm [11"-15.3"]
Footplate Angle	Plantarflexion 10° Dorsiflexion 10°	Plantarflexion 10° Dorsiflexion 10°
Max. Knee Contracture Accommodated	25°	25°
Tray Height Range (from top of hip pad)	0mm-150mm [0"-5.9"]	0mm-150mm [0"-5.9"]
Tray Angle Adjustment	45° - Prone 45° - Supine	45° - Prone 45° - Supine

LECKEY Horizon

Robust, durable power stander
for users aged 4 years to adult.



Product Features

- Accommodation of leg length discrepancy, scoliosis, lumbar lordosis, kyphosis, knee flexion contractures and asymmetrical head position. The frame's growth potential and the wide range of accessories accommodate changing need and potential asymmetry as the user grows.
- Perfectly sized for hoisting and supine standing. The large landing space for hoisted transfers coupled with support pads which provide a proximal fit to suit each user provides the perfect combination.
- Robust, stable and durable frame ensures ongoing confident product use.
- Individually adjustable chest, hip and knee pads as well as footplates provide unsurpassed levels of adjustability.
- Height adjustable hip laterals and pelvic belt, enabling de-rotation of the hips, alignment of the pelvis and centre of gravity over the base of support.



Available in 3 fabric colours

Technical Details			
	Size 1	Size 2	Size 3
Age	4-10	8-14	14-adult
User Height	1000mm-1400mm [39.3"-55.1"]	1400mm-1800mm [55.1"-70.8"]	1550mm-1800mm [61"-70.8"]
Max. User Weight	50kg [110lbs]	80kg [176lbs]	100kg [220lbs]
Shoulder Height (footplate to top of pad)	820mm-1150mm [32.2"-45.2"]	1000mm-1450mm [39.3"-57"]	1000mm-1450mm [39.3"-57"]
Hip Pad Height (footplate to top of pad)	450mm-750mm [17.7"-29.5"]	750mm-950mm [29.5"-37.4"]	750mm-1150mm [29.5"-45.2"]
Hip Pad Width	175mm-270mm [6.8"-10.6"]	210mm-390mm [8.2"-15.3"]	310mm-470mm [12.2"-18.5"]
Knee Height (footplate to middle of pad)	250mm-450mm [9.8"-17.7"]	400mm-550mm [15.7"-21.6"]	400mm-650mm [15.7"-25.5"]
Width Between Chest Laterals	200mm-300mm [7.8"-11.8"]	240mm-380mm [9.4"-14.9"]	340mm-460mm [13.3"-18.1"]
Footplate Angle Adjustment +/-	10°	10°	10°
Angle Adjustment	0°-90°- Prone 0°-90°- Supine	0°-90°- Prone 0°-90°- Supine	0°-90°- Prone 0°-90°- Supine
Tray Angle Adjustment	0°-40°- Prone +/- 20°/30°- Supine	0°-40°- Prone +/- 20°/30°- Supine	0°-40°- Prone +/- 20°/30°- Supine

LECKEY Pronestander

Robust, adjustable prone standing system for users aged 1-18 years.

Product Features

- Flip-up sandals support users completing a standing transfer directly into the product.
- Long base ensures maximum stability for even the most active users, providing caregivers with peace of mind.
- Unique de-rotation belt and pelvic pad work together to keep the user in a central position within the frame with equal pressure over each foot.
- Individually adjustable chest, hip, knee and sandal supports to meet a wide range of clinical needs.
- Easy to set up and adjust with five incremental secure prone positions from 85 to 45 degrees.



Size 1 only available in Meadow Fabric
Size 2 & 3 only available in Mondeo Fabric

Technical Details			
	Size 1	Size 2	Size 3
Age	1-5	4-10	9-18
User Standing Height	800mm-1151mm [31.4"-45.3"]	1100mm-1448mm [43.3"-57"]	1400mm-1849mm [55.1"-72.7"]
Max. User Weight	60kg [132lbs]	80kg [176lbs]	90kg [198lbs]
Chest Pad Height	470mm-680mm [18.5"-26.7"]	680mm-975mm [26.7"-38.3"]	975mm-1350mm [38.3"-53.1"]
Prone Angles	85°, 75°, 65°, 55°, 45°		
Width between Chest Laterals	152mm-254mm [5.9"-10"]	203mm-305mm [7.9"-12"]	254mm-356mm [10"-14"]
Tray Sizes	540mmx400mm [21.2"x15.7"]	590mmx425mm [23.2"x16.7"]	644mmx470mm [25.3"x18.5"]

LECKEY Freestander

Robust, upright standing system for users aged 1-18 years.

Product Features

- Flexibility to alter the frame to meet individual user needs e.g. pelvic support can be removed for users with trunk control, conversely, the tray height and angle can be increased where more support is required. A simple, active set-up is also possible for more able users e.g. hemiplegic or GMFCS II/III.
- Straightforward, intuitive adjustment as all supports run off two vertical uprights.
- Pelvic belt achieves a unique position which enables de-rotation of the hips, alignment of the pelvis and centre of gravity over the base of support as well as accommodation or correction of flexion at the hips and knees.



Size 1 only available in Meadow Fabric. Size 2 & 3 only available in Mondeo Fabric.

Technical Details			
	Size 1	Size 2	Size 3
Age	1-5	4-10	9-18
User Standing Height	800mm-1151mm [31.4"-45.3"]	1100mm-1448mm [43.3"-57"]	1400mm-1849mm [55.1"-72.7"]
Max. User Weight	60kg [132lbs]	80kg [176lbs]	90kg [198lbs]
Chest Pad Height	470mm-680mm [18.5"-26.7"]	680mm-975mm [26.7"-38.3"]	975mm-1350mm [38.3"-53.1"]
Width between Chest Laterals	152mm-254mm [5.9"-10"]	203mm-305mm [7.9"-12"]	254mm-356mm [10"-14"]
Tray Sizes	540mmx400mm [21.2"x15.7"]	590mmx425mm [23.2"x16.7"]	644mmx470mm [25.3"x18.5"]

LECKEY Todstander

Simple and easy to use, upright standing system for kids aged 1-5 years.

Product Features

- Can be used as a stander or by therapists for standing assessments. Promotes active standing, movement and weight shift.
- Compact and easily transported to different environments e.g. from therapy centres to the home.
- Provides a gentle introduction to rehab equipment via a simple, easy-to-use, child-friendly design.

Technical Details	
	One Size Only
Age	1-5
Chest Band Height	470mm- 640mm [18.5"-25.1"]
Base Depth	450mm [17.7"]
Base Width	450mm [17.7"]
Max. User Weight	20kg [44lbs]



Available in 1 fabric colour combination

LECKEY MyWay+

MyWay+ is an open frame walker which cleverly supports children in an upright position, improving stepping and maximising opportunities for exploration, interaction and participation.

A bespoke frame combined with a wraparound hoistable harness and a wide range of accessories, ensures the MyWay+ can be tailored to meet every need and ability.



Product Features

- The unique corset-like harness gives core trunk stability for optimal stepping pattern and active movement. The MyWay+ harness provides the support required for users to develop head and limb control. Its hug-like effect increases sensory feedback, dampens involuntary movement and increases confidence and the chance of succeeding.
- MyWay+'s open-front, hands-free frame maximises opportunities for exploration and participation. The open-design also allows therapists to get close to the child and work with their lower limbs whilst they are in MyWay+, facilitating guided movement.
- The hoistable harness makes transfers easier- direct transfers eliminate the need for additional hoisting slings and make the most out of therapy time. Hoisting in the MyWay+ harness also provides superior comfort to the user as it does not slip into the armpits.
- Castors – Indoor and Outdoor (Robust, anti-roll back, direction lock [including *sideways lock], varied resistance) – to suit a child's individual needs.
- Range of optional accessories available to meet individual requirements, including headrest, ankle guides and saddle.




Each harness size is a different color

Technical Details			
	Size 1 Frame	Size 2 Frame	Size 3 Frame
Age (approx.)	1-5	4-11	10-16
User Shoulder Height (to floor)	590-870mm	870-1190mm	1070-1490mm
User Saddle Height (inside leg)	240-410mm	400-610mm	550-810mm
Max. User Weight	25kg	50kg	80kg
Frame Weight	14kg (Size 1 frame and upper support unit do not separate)	7.5kg	11kg
Upper Support Unit Weight		10kg	11kg
Frame Width (mm)	570mm	670mm	720mm
Frame Length (mm)	840mm	930mm	1180mm

Harness						
	Size 1	Size 2	Size 3	Size 4	Size 5	Size 6
Age (approx.)	1-4	2-6	5-9	8-12	12-15	Up to 16
Waist at Navel Level (cm)	46-56	50-60	55-65	60-70	70-80	80-90
Colour	Purple	Green	Blue	Red	Orange	Grey

LECKEY MyWay Pedal

MyWay Pedal is a revolutionary elliptical trainer which helps improve strength and functional ability. Attaching to the size 2 and 3 MyWay frame, MyWay Pedal is ideal for giving on-the-go kids the opportunity to be active, have fun and join in with their friends.

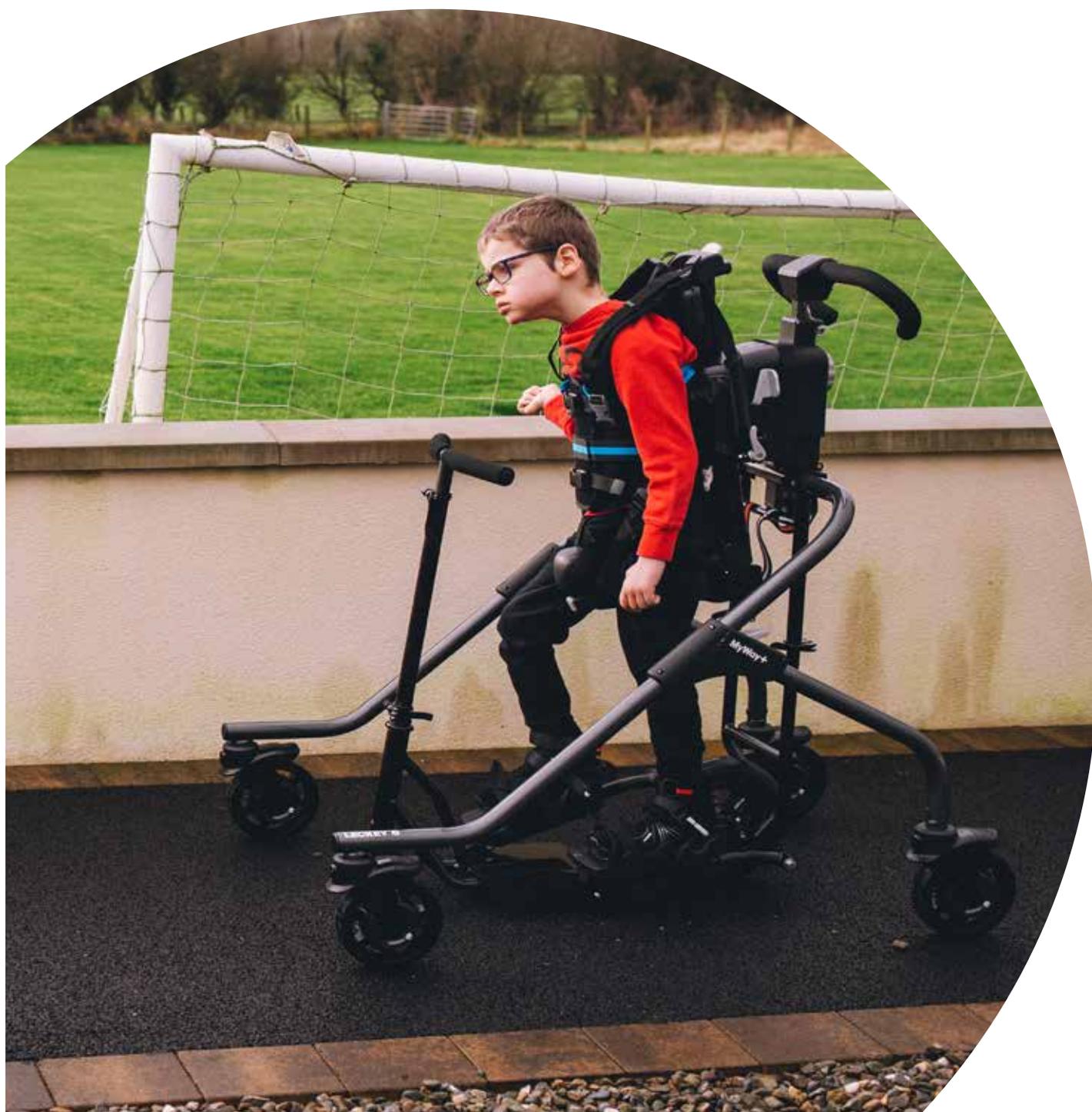

Product Features

- Targets and strengthens anti-gravity muscles to help reduce the risk of hip and knee contractures and improve function. Pedalling strengthens anti-gravity muscles (calf, quads and hamstring)- into extension unlike a trike which encourages a flexed position. MyWay Pedal ensures the quads are worked exactly where strength is often lacking, at the end of range. Working the end of range will reduce the risk of hip and knee contractures which may limit standing and walking ability.
- Facilitates elliptical action more like stepping to help support gait therapy goals. MyWay Pedal provides a sensori-motor experience similar to gait. Unlike the circular motion of cycling, the elliptical pedalling of MyWay Pedal gives an extended stretch of the lower limbs as it involves a longer forwards/backwards movement rather than an up/down movement. It helps improve loading around the hip to aid joint development.
- Use in conjunction with the Leckey Connect App to chart performance and motivate the child.



The app provides a valuable way of measuring performance and progress on a child's speed, distance and duration of elliptical pedalling. It is available on Android Smartphone, tablet and the Leckey online portal.

- Attaches to the size 2 and 3 MyWay frame and can be fitted with the user in or out of the frame.
- Two sandal sizes available.
- Height adjustable and removable handlebar.



LECKEY



Other Products

LECKEY Advanced Bath Chair

The Leckey Advanced Bath Chair provides postural support while ensuring the individual is safe and secure during bathing.

Product Features

- Provides the user with appropriate levels of postural support ensuring they are safe and secure at all times during bathing.
- Lightweight and robust frame.
- Simple and easy-to-use adjustments mean that the product can be set up in a variety of positions in seconds.
- Height adjustable legs.
- Breathable cover, removable and machine washable.
- Back angle adjustment from 0-90 degrees.
- Adjustable hip belt.
- Chest belt with integrated lateral supports.
- Easy to use hand adjustments.
- Folds flat for storage.



Technical Details	Size 1	Size 2	Size 3	Size 4
Age	1-5	4-9	8-14	12-18
Min User Height	750mm [30"]	950mm [37"]	1150mm [45"]	1350mm [53"]
Max User Height	1050mm [41"]	1250mm [49"]	1450mm [57"]	1650mm [65"]
Max User Weight	72kg [159lbs]	72kg [159lbs]	72kg [159lbs]	72kg [159lbs]



Available in 2 fabric colours

LECKEY Shower Trolley

The Shower Trolley complements the Advanced Bath Chair which can be converted within seconds into a shower chair.

Product Features

- Converts the Advance Bath Chair to a versatile shower chair in seconds.
- Manufactured from medical grade stainless steel.
- Robust frame and high quality lockable castors.
- Can be easily and safely manoeuvred over wet and slippery surfaces.
- Can be adjusted to fit all four sizes of the Leckey Advance Bath Chair.



Technical Details	
One Size Only	
Age	1-18
Length	900mm-1150mm [36"-45"]
Width	600mm [24"]
Seat Height from Floor	650mm-700mm [26"-28"]

LECKEY Sleepform

Strong, durable, practical and adjustable sleeping system for kids aged 0-18 years.

Product Features

- Close proximal fit cocoons the child and makes them feel safe and secure. Adapts quickly and easily to changing body shape.
- Compact and portable for visits to family, going on holiday etc. Suitable for use on a variety of surfaces.
- Supports 24-hour postural management as can be used during the day to support changes of position.
- Can be used in cots and single beds and suits all sleeping positions.
- Air flow mattress.
- Grid sheet to record best sleeping positioning.
- Temperature control sheet available.




Technical Details				
	Size 1	Size 2	Size 3	Size 4
Age	0-1	1-5	1-5	5-18
Suitable for	Cot	Cot	Single Bed	Single Bed

Supporting Children





Phone: 0800 31 61 81
Email: helpis@alliedmedical.co.nz
www.alliedmedical.co.nz

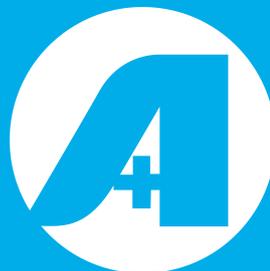

**ALLIED
MEDICAL**
LIMITED