
allard
Support for better life!

Instructions for Use

Vission Finger Splint POSI

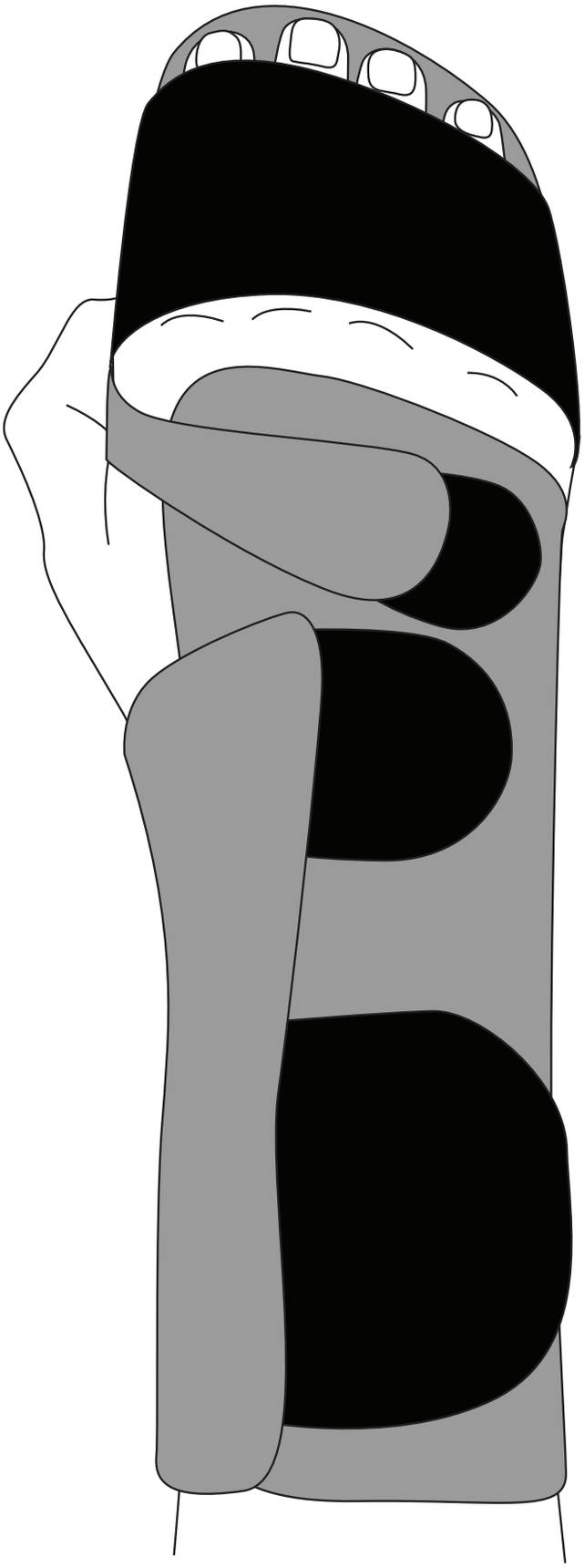

EN Vission Finger Splint POSI

46019 Vission Finger Splint POSI
Carefully read these instructions before use!

Product Description
A textile volar splint with an adjustable aluminum stay covering fingers II-V, palm and wrist, positioning the joints in Position Of Safe Immobilization (POSI).

Intended Use
For temporary immobilization of the wrist and fingers following trauma or surgery and/ or preventing contractures of the finger joints.
Single patient - Multiple use.

Indications
Trauma, Post-Op, fractures of metacarpals and fingers. Dupuytren's contracture.

Contraindications
No known

Limitations
Not suitable for persons with neurological disorders that may include impaired sensitivity, spasticity, and/or contractures.

Adjustments
Fitting and modifications of this product may only be performed by qualified healthcare professionals.

Warning
The patient should check their skin daily. If there is any sign of skin irritation, significant swelling or unexpected severe pain, the orthosis should be removed immediately.

If the patient has been instructed not to remove the orthosis during treatment, healthcare professionals should be consulted first.

Material specifications
Polyamide, Lycra, Polyurethane, Polyester, Aluminium. The product is not flame retardant.

Washing instruction
Close hook and loop and remove stay prior to laundry. Use a washing bag. Note the position of the stay to ensure correct position after washing.

Machine wash at 40° C, gentle cycle. Dry flat. Do not bleach, tumble dry, iron or dry clean.

Questions
For any questions or issues, please contact your healthcare professional/fitter.

Disposal
Separate the metal stay from the textile and dispose each part in accordance with local legislation.

Information
Any serious incident that may occur in relation to the device should be reported to the manufacturer and the competent authority of the Member State in which the user and / or patient is established.

Instructions


Check that the product is complete and without damage!

1. **Positioning:** Insert the hand into the orthosis and secure the straps, starting with the one across the back of the hand to stabilize the splint.


2. **Adjustments of finger strap:** If needed, adjust the position and the length of the finger strap.


3. **POSI**: The stay is pre-shaped in a "Position Of Safe Immobilization" (POSI). Adjustments to the stay should only be performed by qualified healthcare professionals.

1

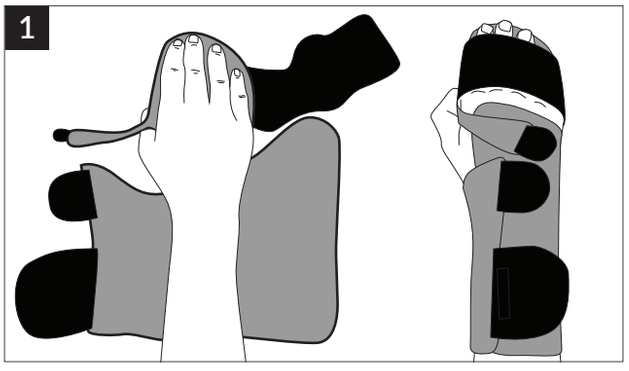


2

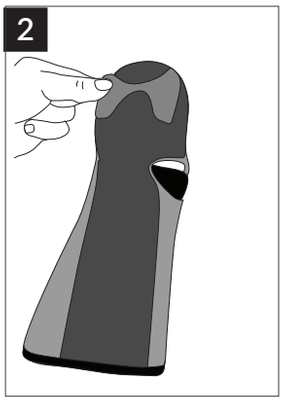

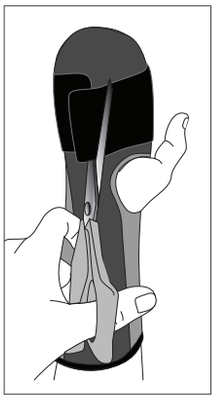

3

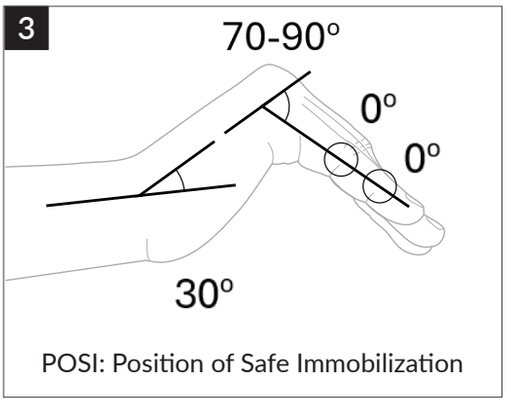

POSI: Position of Safe Immobilization

Art.no.	Size	Left/Right	Wrist circumference
460191011	S	Left	13-16 cm
460191012	M	Left	16-19 cm
460191013	L	Left	19-22 cm
460191014	XL	Left	22-25 cm
460192011	S	Right	13-16 cm
460192012	M	Right	16-19 cm
460192013	L	Right	19-22 cm
460192014	XL	Right	22-25 cm





ALLIED MEDICAL