

# BACKUP



REHASENSE

## SHOWER CHAIRS



BACKUP - a powder-coated, height adjustable aluminium shower chair with fixed handles and an optional backrest.

## • HEIGHT ADJUSTABLE

Seat height 44-64 cm, can be adjusted to 8 different positions with 2,5 cm steps by the use of snaplock clips, so it is suitable also for tall users. Seat width and depth: 40×40 cm.

## • REMOVABLE BACKREST

Removable and length adjustable backrest is an accessory on the fixed frame showerchairs and not installed as standard.

The backrest can be adjusted to 2 different positions in increments of 1 cm. Chairs can be used with or without a backrest.


## • DURABLE & STRONG

The construction has been tested and approved on users with a maximum weight of 150 kg.



The removable seat and backrest is made from skin-friendly, anti-slip polyurethane with a textured surface that ensures optimal drainage of water and soap residue.

The blow-moulded HDPE material warms quickly to body temperature under the shower.

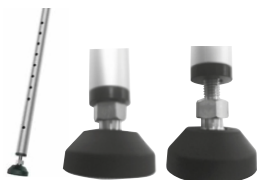

With one fine tuning leg for final, precise adjustment and levelling of the chair.



All legs with 8×2,5 cm increments to adjust seat height between 44-64 cm. Rubber ferrules prevent the chair from slipping on a slippery surface.

### SPECIFICATIONS

Max user weight:  
150 kg

Colours:  
Frames: White or grey  
Seat & back: Grey

Recommended user height:  
130-210 cm

Width×Length×Height:  
56×53×60-80 cm (without backrest)  
56×59-60×80-100 cm (with backrest)

Width between armrests:  
42 cm

Seat height:  
44-64 cm (every 2,5 cm)

Total weight:  
4,8 kg (incl. backrest at 2,6 kg)

Material:  
Non-corrosive, powder coated  
aluminium & PU (polyurethane)