

ICON 120

THE COMFORT WHEELCHAIR



ICON 120

COMFORT WHEELCHAIR

MODULAR THINKING

As part of the **ICON Mobility System**, utilising modular design and interchangeable parts across the ICON family, the **ICON 120** can cater for many different user types and functional needs. All day comfort is ensured with the reclining back and tilt in space seat when combined with pressure distributing seat and back cushions. The arm rests and head rests are covered in soft and pliant PU fabrics to minimise pressure concentrations.

TECHNICAL INFORMATION

Safe weight limit:	180 kg
Safe weight limit for use as seat in car:	134 kg
Weight:	39 kg
Transportation weight:	27 kg
Colors:	white or sober grey



Frame sizes for best fit: 32, 36, 41, 46, 51, 56 and 61 cm.



The chair can be set up as Self propelled with 20", 22" and 24" rear wheels.



For the Attendant Propelled set up, the rear wheel diameter is 16".



4-point tie down system for transport and securing the wheelchair in a vehicle.



The wheel base can be extended at the front by 5 and 10 cm.




Brake options include attendant drum brakes and knee brakes with short or long lever.



Height adjustable armrest in a range of 27 – 40 cm.



The wheelchair is equipped standard with height adjustable push handles and multi-adjustable head rest.

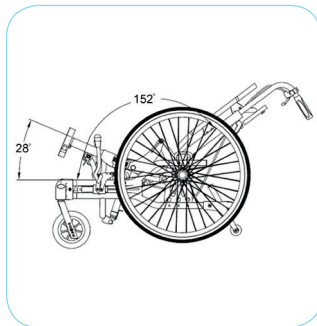

152° total tilt and recline.



Front wheel options include puncture proof PU in 6", 7" and 8" diameters at multiple positions in the front fork.



Lateral supports as an option.



The wheelchair is equipped standard with length adjustable elevating legrests.