

Let's Shower

Shower chair

Black TCLSH101
White TCLSH201



Let's Shower is a new shower stool in a modern design intended for use in the shower or the bathroom. The stool features a beautiful and stable construction for safe and secure use. The legs are easy to adjust to desired height with one easy motion. The stool is easy to fit in any space thanks to its small shape.

Let's Shower is adjustable from 45 to 53 cm in height. The height can be adjusted by pressing the buttons on the legs and then moving the lower tubes to the preferred height, the buttons will then pop back into place and lock the legs.

Important: Always adjust all the legs to the same height to avoid falling accidents.



White. Art. No. LSH201


Technical specifications:

Tested according to ISO 12182:2012

Max. user weight 150 kg
Adjustable height 45 – 53 cm
Weight 2 kg



Material:
Seat ABS
Tube frame Aluminum
Feet TPE

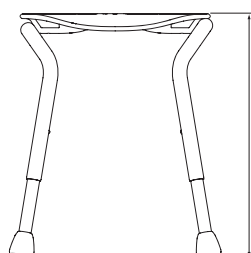

45 - 53 cm

