




CE MD MEDICAL DEVICE - CLASS I

mewa

SPECIAL NEEDS STROLLER FOR CHILDREN

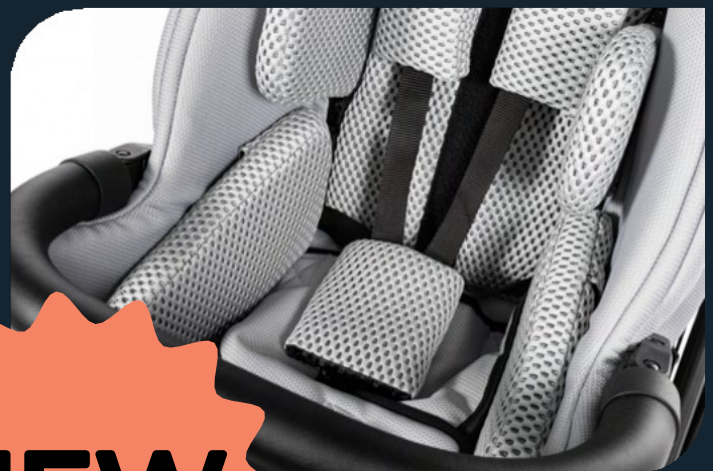

NEW

GENERAL INFORMATION

The **MEWA special needs stroller** is designed for the youngest users who are unable to walk or sit in a **MEWA** as a result of injuries, diseases or motor dysfunctions. The stroller is suitable for children at a height of **90 – 125 cm** and a maximum weight of **40 kg**.

Thanks to the backrest tilt and recline feature as well as adjustable torso and head pads you may easily adjust **MEWA** to the required position and ensure the user with optimal protection and stability.

The shock-absorbing suspension on all four, large, gel-filled wheels together with the tilt-foot facilitate travel on uneven terrain. Modern and stylish design is an additional advantage of this stroller.

INDICATIONS

- paralysis and paresis of lower limbs various origins
- paralysis and paresis of the trunk
- weakness of postural muscles, responsible for maintenance correct body position
- cerebral palsy
- myelomeningocele
- muscular dystrophies

CONTRAINDICATIONS

In case of any contraindications to sitting position &/or passive verticalization stroller adjustments should be consulted with physiotherapist or physician-in-charge.



LIGHTWEIGHT CONSTRUCTION

The lightweight design of the stroller is an ideal solution for parents who value mobility and comfort. The stroller is easy to carry and maneuver, which makes walks and trips less burdensome. At the same time, it also makes it easier to fold and store the stroller.



SHOCK ABSORPTION

Shock absorption in strollers is crucial for comfort and safety. It effectively dampens shocks on uneven terrain, which allows the child to travel in peace and without any discomfort. At the same time, it makes it easier to drive the stroller, allowing smooth maneuvering, which reduces the burden on the parent.



SHOPPING BASKET

The shopping basket with a capacity of up to 3 kg, mounted on the **Mewa** stroller, is a convenient solution for parents. It allows you to easily transport necessary accessories and small items, without the need to carry additional bags.



UV 50 TECHNOLOGY

UV technology in fabrics protects against fading, so colours remain intense despite exposure to the sun.



NEW UPHOLSTERY MATERIAL

The material used in the construction of the stroller is friendly to both people and animals. Made of safe, non-toxic raw materials, it meets the highest quality and ecological standards. It is gentle on the child's skin, and at the same time durable and resistant to damage.



BREATHABLE MATERIAL

The interior of the stroller and part of accessories are made of breathable material that provides excellent air circulation. Thanks to this, the child feels comfortable even during long walks, and the skin remains dry and free from irritation. The material is not only breathable, but also soft and delicate, which additionally increases the comfort of use.



DURABLE POLYURETHANE TIRES

Extremely durable and flexible tires, made entirely of polyurethane, so there is no need to pump them up. This design significantly reduces the risk of accidental damage, ensuring reliability and long service life.



REVERSIBLE SEAT

The design of the stroller allows you to easily rotate the seat opposite to the driving direction. Thanks to this feature you can keep an eye on your child all the time.



CANOPY FOLDED




CANOPY UNFOLDED




CANOPY VENTILATION

The new design of the canopy for the stroller combines elegance and functionality. Thanks to the innovative system of extending the visor, the canopy provides even better protection from the sun, adapting to changing weather conditions. Additionally, the built-in ventilation guarantees optimal air flow, which is crucial on hot days, protecting the child from overheating.

OTHER UPHOLSTERY COLORS

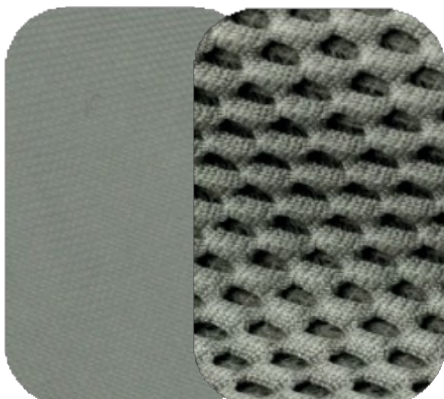

GREEN

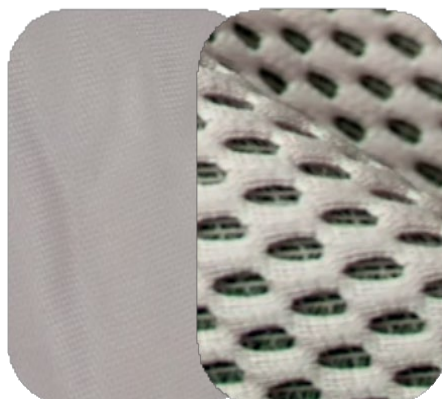

BEIGE



BLACK

EQUIPMENT



SET MINIMIZING KIT



HEAD/TORSO
ADJUSTABLE PADS



ABDUCTION BELTS



SAFETY BARRIER



FIVE-POINTS BELTS



ADJUSTABLE
WEDGE



THERAPY TRAY



CANOPY WITH
WINDOW



WINTER GLOVES




TRAVEL BAG



SHOPPING BASKET



SUN UMBRELLA

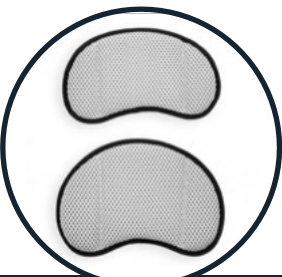

CUSHION PILLOWS



SLEEPING BAG



MOSQUITO NET



RAIN COVER

TECHNICAL PARAMETERS OF THE STROLLER

seat width: 25*/32 cm
effective seat depth: 24*/28 cm
backrest height: 62 cm
backrest width: 20*/30 cm
distance between footplate and seat: 17 - 32 cm, 5 positions

*values are given with pads

total weight: 20,5 - 25 kg
stroller length: 97 - 132 cm
stroller height: 125 - 128 cm
stroller width: 62,5 cm



IMPORTANT NOTE! The choice of a medical device appropriate for a specific disease entity and condition of the patient must be made with the help of a specialist doctor or physiotherapist. Before use, read the operating manual and the label or consult your doctor, as this medical device may not be suitable for you.