


mouse

SPECIAL NEEDS STROLLER FOR CHILDREN

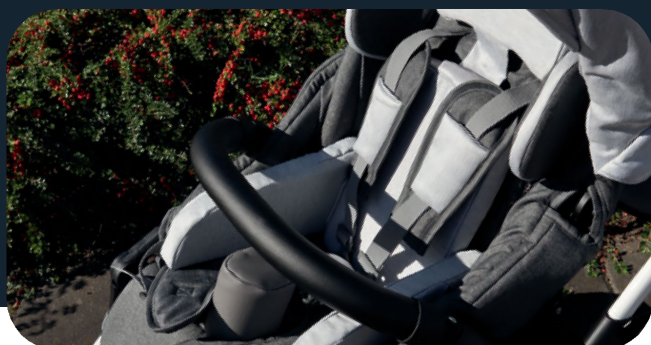

CE **MD** MEDICAL DEVICE - CLASS I

GENERAL INFORMATION

The **MOUSE** special needs stroller is designed for the youngest users who are unable to walk or sit in a correct position as a result of injuries, diseases or motor dysfunctions. The stroller is suitable for children at a height of **80 – 110 cm** and a maximum weight of **30 kg**.

Thanks to the backrest tilt and recline feature as well as adjustable torso and head pads you may easily adjust **MOUSE** to the required position and ensure the user with optimal protection and stability.

The shock-absorbing suspension on all four, large, gel-filled wheels together with the tilt-foot facilitate travel on uneven terrain. Modern and stylish design is an additional advantage of this stroller.

INDICATIONS

Paralysis and paresis of lower limbs of various origins, paralysis and paresis of the trunk, weakening of postural muscles responsible for maintaining the correct body position, cerebral palsy, myelomeningocele, muscular dystrophies.

CONTRINDICATIONS

In case of any contraindications to sitting position &/or passive verticalization stroller adjustments should be consulted with physiotherapist or physician-in-charge.

FEATURES / SAFETY AND COMFORT

TILT IN SPACE SEAT

simple mechanism



8° - 23°

4 positions

RECLINING BACKREST

wide range of adjustments



16° - 88°

4 positions

ADJUSTABLE AND ANGLE TILT FOOTREST

simple and quick adjustment



84° - 173°

6 positions

SEAT MINIMIZING KIT

reducing insert



The **seat minimizing kit** in **Mouse** is a comfortable, one-piece insert that allows you to reduce both width and depth of the seat.

STABILIZATION SYSTEM

adjustable pads



The **head and torso support** consists of two pairs of independent pads that are height, width and angle adjustable. Its wide range of adjustments and simplicity gives you a fast and effective adaptation to the user's individual needs.

LYING POSITION

wide range of backrest recline



With a fully adjustable seat and footrest you may quickly alter the stroller into cozy bed, perfect for a short rest or a nap.

REVERSIBLE SEAT

arrange your travel easily

FOLDING LIGHTWEIGHT FRAME

ease of storage, transportation
and dismantling

DRIVING SUPPORT

comfort for the bumpiest journeys



The design of the stroller allows you to easily rotate the seat opposite to the driving direction. Thanks to this feature you can keep an eye on your child all the time.



The construction of main elements allows you to dismantle **Mouse** into 3 separate units (frame, seat and wheels).



The stroller is equipped with **a shock-absorbing suspension** on all four wheels and **an adjustable rear suspension** in the rear wheels with option to lock, **central brake** and **a tilt-foot**.

FUNCTIONAL ACCESSORIES IN PRACTICE

**1**

Canopy with window



Set minimizing kit

**2**

Torso/head pads

**4**

Adjustable footrest

**5**

Swivel front castor wheels with option to lock

**6**

Quick release mount on all 4 wheels



EQUIPMENT



SET MINIMIZING KIT

HEAD/TORSO
ADJUSTABLE PADS

ABDUCTION BELTS



SAFETY BARRIER



FIVE-POINTS BELTS

ADJUSTABLE
WEDGE

THERAPY TRAY

CANOPY WITH
WINDOW

WINTER GLOVES



TRAVEL BAG



SHOPPING BASKET



SUN UMBRELLA



RAIN COVER



SLEEPING BAG



MOSQUITO NET

TECHNICAL PARAMETERS OF THE STROLLER

seat width: 24*/32 cm
 effective seat depth: 22*/26 cm
 backrest height: 52 cm
 backrest width: 30 cm
 distance between footplate and seat: 14,5 - 23 cm, 4 positions
 *values are given with inserts

total weight: 17,2 - 23,2 kg
 stroller length: 95 - 124 cm
 stroller height: 122 - 134 cm
 stroller width: 60 cm



IMPORTANT NOTE! The choice of a medical device appropriate for a specific disease entity and condition of the patient must be made with the help of a specialist doctor or physiotherapist. Before use, read the operating manual and the label or consult your doctor, as this medical device may not be suitable for you.