

PIPER COMFORT



This rehab buggy manufactured by Patron® is a state of the art mobility solution for children.

With premium design at affordable prices, this buggy is great for light outdoor activity, whilst being fully crash tested for transportation.

This robust product corresponds to the newest MDR + EC standards for hygiene, technical and safety requirements.


FUNCTIONS INCLUDED

- > Extended backrest recline
- > Extended backrest height
- > Large adjustable footrest
- > Super compact foldable design
- > Load capacity up to 50kg

PARTS INCLUDED

- > Piper Comfort Buggy
- > Premium **NEW** Upholstery
- > 5 Point Immobilising Buggy Harness
- > Adjustable One Piece Footplate
- > Shopping Basket
- > Front Wheels with Direction Lock
- > Brakes with transport lock
- > Crash Tested Tie Down Kit
- > Detachable Wheels for Easy Folding



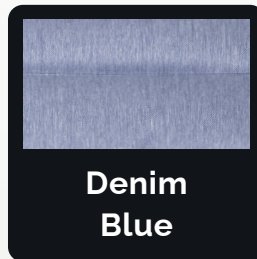

SIMPLE SEAT REDUCTION

The Piper Comfort has optional extras of a seat inlay, set of width reducers and a soft footrest raiser to reduce the size of the buggy for smaller children.

This provides great growth in the system as the reducers can be easily removed as the child grows!





PIPER DIMENSIONS

	Standard	With Inserts
Backrest Height	680 mm	680 mm
Seat Depth	290 mm	240 mm
Seat Width	260 to 340	230 to 270
Seat to Footrest	240 to 370	200 to 230
Backrest Angle	90° to 125°	90° to 125°
Seat Angle	25°	25°
Total Length	1080 mm	1080 mm
Total Width	590 mm	590 mm
Total Height	1000 mm	1000 mm
Max User Weight	50.0 kg	50.0 kg

FOLDED DIMENSIONS

	Standard	With Inserts
Buggy Weight	14.5 kg	16.0 kg
Folded Length	1020 mm	1020 mm
Folded Width	330 mm	330 mm
Folded Height	280 mm	280 mm



ENCLOSED BACKREST RECLINE

The Piper Comfort features a unique backrest recline system, offering 35° of movement whilst keeping the buggy user completely enclosed within the frame.

This recline system is great for users with problems relating to sensory overload or behavioural issues, as the buggy cocoons their body, making them feel safe and secure whilst out on the move.

The recline also brings up the footplate, allowing the user to almost lie flat.



A RANGE OF ACCESSORIES

