


## Change of the anti-tipper on Flexx series

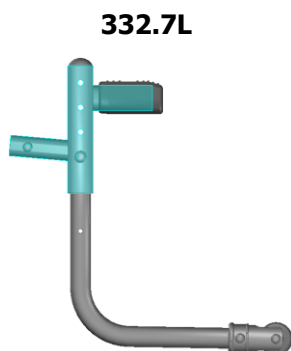

Previous Version	Updated version
<p><b>109mm</b></p>	<p><b>129mm</b></p>

- ✔ In response to the anti-tipper spring accidentally detaching when a user was operating it by foot, Karma has changed the length of the outer tube (note in red) from 109mm to 129mm. Note that this change will NOT affect the total length of the wheelchair.
- ✔ The updated version is backward compatible with the existing one.
- ✔ The updated version was implemented by the batch of serial no. TM22-074041.
- ✔ The updated version was implemented on Flexx series.

### Order Information

Article number	151200000152	KM8022 LeftRetractile Anti-Tipper Assy 332.7L
	151200000100	KM8022 RightRetractile Anti-Tipper Assy 332.7L
	151200000099	KM8022 LeftRetractile Anti-Tipper Assy 118L
	151200000153	KM8022 RightRetractile Anti-Tipper Assy 118L