

Explore your possibilities


- ▶ Stays with you during activities
- ▶ Skin protection outside the wheelchair
- ▶ Extremely light to support optimal mobility

Ergonomically shaped back part for seamless fit to the shape of your lower back



Easily adjustable straps to provide comfort and to follow the body

Vicair SmartCells for skin-protection

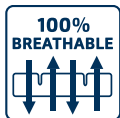

The Vicair AllRounder O2 is a cushion that is strapped on to your body. It provides a moderate level of skin protection for individuals who are at risk of developing pressure injuries.

This occasional use activity cushion is used to provide comfort and prevent pressure injuries in situations outside the wheelchair. The Vicair AllRounder O2 consists of 5 individual compartments, designed to contour under loading. Each compartment is designed to influence a specific area of the pelvis. The standard filling of the AllRounder O2 is ready to use out of the box and fits most users who: play/sit on the ground, sail, handbike, go up and down stairs, doing chores, etc.

Product performance

Skin protection / Anti decubitus

Pressure redistribution	◀◀◀◀◀◀
Shear force reduction	◀◀◀◀◀◀
Microclimate control	◀◀◀◀◀◀
Impact & vibration damping	◀◀◀◀◀◀

Positioning

Lateral stability	N/A
Frontal stability	N/A
Sliding reduction	N/A
Allows manual adjustments	N/A
Increasing endurance	N/A

Maintenance & Safety

Durability	◀◀◀◀◀◀
Reliability	◀◀◀◀◀◀
Functionality when misused	◀◀◀◀◀◀
Easy to maintain	◀◀◀◀◀◀
Level of hygiene	◀◀◀◀◀◀

Functionality & application

Comfort	◀◀◀◀◀◀
Ultra lightweight	◀◀◀◀◀◀
Allows easy transfer	◀◀◀◀◀◀
Available standard sized	◀◀◀◀◀◀
Maximum user weight	◀◀◀◀◀◀

Technical specifications

Thickness	N/A
Average weight	800 grams
Max. user weight	100 kg
Machine washable	60°C max.

Cushion sizes

Small	< 38 cm
Medium	28 - 42,5 cm
Large	43 - 48 cm

Determining the correct hip width

Measure the width of your hips to check which size matches the size of your body. When seated, hold a tape measure horizontally above the widest point of your upper legs, closely against the front of your hips. Do not follow your leg contours, keep the tape measure straight. Measure from A to B (or B to A). Choose the larger size if the measured size is in between two sizes.

figure 1

