

Offloading pressure in the rectal area

- ▶ For completely offloading of the rectal area
- ▶ Adjustable non-contact area
- ▶ Stable positioning



PRESSURE INJURY TREATMENT PLAN

SOLID SHAPE

READY FOR USE

100% BREATHABLE

MACHINE WASHABLE

The Vicair Centre Relief O2 consists of 5 individual compartments, designed to contour under loading. Each compartment is designed to influence a specific area of the pelvis, striking a unique balance between preserving skin and underlying tissue at the ischials and trochanters, positioning the pelvis and supporting the femurs. The unique feature is the circular compartment towards the rear of the cushion. Designed especially for (partial) offloading of the coccyx/rectal area. For specific needs the filling of each of the 5 compartments can be easily adjusted, for individuals who: (1) require (partial) offloading of the rectal area, (2) have coccyx problems, (3) have localized deep pressure injuries in the designated area.

Product performance

| Skin protection / Anti decubitus | |
|----------------------------------|-------|
| Pressure redistribution | ◀◀◀◀◀ |
| Shear force reduction | ◀◀◀◀◀ |
| Microclimate control | ◀◀◀◀◀ |
| Impact & vibration damping | ◀◀◀◀◀ |

| Positioning | |
|---------------------------|-------|
| Lateral stability | ◀◀◀◀◀ |
| Frontal stability | ◀◀◀◀◀ |
| Sliding reduction | ◀◀◀◀◀ |
| Allows manual adjustments | ◀◀◀◀◀ |
| Increasing endurance | ◀◀◀◀◀ |

| Maintenance & Safety | |
|----------------------------|-------|
| Durability | ◀◀◀◀◀ |
| Reliability | ◀◀◀◀◀ |
| Functionality when misused | ◀◀◀◀◀ |
| Easy to maintain | ◀◀◀◀◀ |
| Level of hygiene | ◀◀◀◀◀ |


| Functionality & application | |
|-----------------------------|-------|
| Comfort | ◀◀◀◀◀ |
| Ultra lightweight | ◀◀◀◀◀ |
| Allows easy transfer | ◀◀◀◀◀ |
| Available standard sized | ◀◀◀◀◀ |
| Maximum user weight | ◀◀◀◀◀ |

Technical specifications

| | |
|------------------|-----------|
| Thickness | 10 cm |
| Average weight | 750 grams |
| Max. user weight | 250 kg |
| Machine washable | 60°C max. |

Cushion sizes

- Standard
- On request

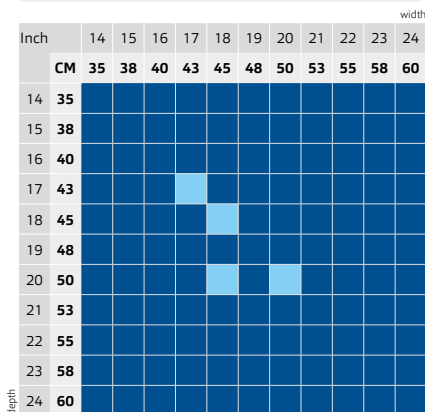

Available covers

Vicair O2 Top Cover, Comfair Cover and Incotec Cover.