

Move more.  
Live more.

Compact, powerful and versatile, the WHILL Model C2 accompanies you on all your adventures, indoors and out.



ALLIED  
MEDICAL  
LIMITED







## Intuitive and manoeuvrable Convenient, even in tight spaces



**Retractable arm rest**  
Lift up the arm rest to easily sit down in the chair from the side. Stand on the foot rest for added stability.



**76 cm turning radius**  
The WHILL is easy to manoeuvre and can navigate tight spaces, making it ideal for use in small indoor areas.



**Modern, intuitive controller**  
Turn the power on and control the speed with one hand. The WHILL will automatically brake when you release the controller, even on an incline.  
The controller can be positioned on either the left or right side of the device.







## Powerful, superior handling

The WHILL Model C2 takes you places where other wheelchairs can't.

Thanks to its two powerful motors and innovative wheels, the WHILL Model C2 can be used on terrain that a standard wheelchair can't handle. Holes in the road, uneven terrain, grass and gravel will no longer be an issue.



Effortlessly negotiate obstacles up to 5 cm high.



Rear wheel suspension absorbs the impact when driving over bumpy roads or obstacles.



The WHILL Model C2 is designed to prevent drifting. Drive straight on a tilted path.





## Compact and easy to transport Take it with you wherever you go!

The WHILL Model C2 is lightweight, can be quickly disassembled and fits into the boot of any car. It can also accompany you on your flight, in the baggage hold.



### How to disassemble the WHILL

Find out how to disassemble the device in four steps by clicking this link and taking a photo of the QR Code:







5 hour  
charge

Range  
18 km

Weight  
2.7 kg

### Removable Lithium-ion battery

Guarantees power and autonomy all day long.

Developed in partnership with LG, the battery charges in five hours, providing the WHILL Model C2 with a range of 18 km, and enough power to last all day. Optional batteries and chargers are available for purchase separately.

The battery can be removed and charged by plugging in to a regular household outlet.



### Control your WHILL with the app\*



- > Control the WHILL remotely.
- > Set the maximum speed, as well as sensitivity of the controls and acceleration.
- > Monitor the battery in real time.
- > Track your activity: travel distance and driving time.

\*WHILL Model C2/F Series





# Choose the colour that suits you

The arm covers are available in ten carefully selected colours. Easily change colour to suit your mood or personal style!



White



Blue



Black



Light Green



Red



Pink



Light Blue



Gold



Grey



Purple



## Testimonials

« Above all, my WHILL is my best friend, allowing me to regaining my freedom of movement! » Olivier



« I wanted a wheelchair that doesn't look like a wheelchair. I wanted it to be a real eye-catcher. And there's so much you can do with it, it's incredibly innovative! I'm much more active now that I've got my WHILL. And when people see me and my wheels, they call out so I can have a chat with them. » Nicolas



« After trying out a number of electric wheelchairs, my favourite way of regaining my freedom, I chose the WHILL C2 without a moment's hesitation! With its unrivalled stability on slopes and extremely smooth, responsive ride, it negotiates obstacles with ease and effortlessly tackles potholes. Yes, you can definitely say I've regained my freedom. » Bertrand

« Faced with a disability, I was so grateful I discovered WHILL. It helped me rebuild my life. It looks good, is easy to manoeuvre in small spaces, copes with any terrain and can clear 5 cm obstacles thanks to its innovative wheel technology. It's my personal vehicle, and it's enabled me to get back on the move and enjoy life to the full again, with those who are dear to me. The WHILL is a compact design wheelchair. It evokes dynamism and movement is life! » Valérie




« At first, I thought the remote control was a gimmick, but it has become really useful. I use it to move from my bed to an armchair, to take a shower, to move the WHILL and put it away, and for visits to the doctor, etc. » Pascal

« As I'm quite small, I needed a lightweight, manoeuvrable solution that could also be folded up for easy storage in the car. The WHILL has given me more autonomy and independence, and for me that's really important. » Stéphanie

# Technical specifications

The Model C2 can be easily adjusted from seat height and back support angle to controller position (right or left drive) to provide a comfortable driving experience.

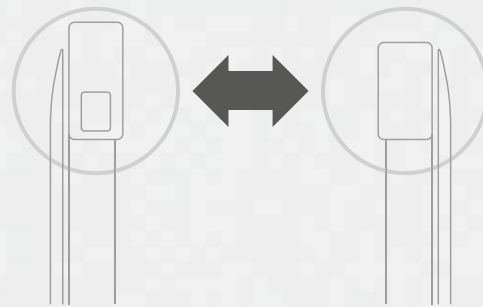

The controller can be positioned on either the left or right side

1. Controller Height .....	670, 690, 710, 730 mm
2. Back Support Angle.....	90, 95, 100°
3. Back Support Height .....	340 - 460 mm
4. Seat Height.....	405, 425, 445, 465, 485 mm
5. Seat Angle .....	5°
6. Seat Depth.....	400 - 450 mm
7. Seat Width.....	400 - 500 mm
8. Foot Support Angle.....	2°
9. Foot Rest Height.....	127 mm

Warning: When using your own cushion or back support, for safety, be sure to use ones that are sufficiently flame retardant. Please note that this brochure is not a complete user manual. Before using the product for the first time, it is essential to refer to the user manual.



Driving Range .....	18 km	Weather Tested .....	IPX4
Maximum Incline .....	10°	Maximum Speed .....	6 km/h
Turning Radius.....	760 mm	Obstacle Clearance.....	5 cm
Device Width .....	554 - 654 mm	Device length.....	985 mm
Weight Fully Assembled .....	51.9 - 52.8 kg	Suspension.....	4 independent wheels
Air Travel Approved .....	Yes	App control.....	Yes, remote control
Charging Time.....	5 hours	Interchangeable arm covers ...	10 colours
Weight Capacity .....	40 cm: 115 kg, 45 cm/50 cm: 136 kg	Battery.....	Lithium-ion (25.3V 10.6Ah)

Note: Driving range of 18 km will vary due to user weight, drive surface, terrain, ambient temperature, speed setting, battery % and its condition.

WHILL Model C2 is classified as Wheelchair Class A.