

karma VIP2

This well optioned tilt-in-space and reclining wheelchair provides the optimal pressure relief and positioning.



ALLIED
MEDICAL
LIMITED



Foldable & lightweight for easy travel

+ FEATURES

Tilt-in-Space

The Tilt-in-Space system allows the wheelchair user to change position throughout the day, which will redistribute and relieve pressure.

Anti-Sliding Reclining System

Sliding is a major concern for reclining chairs. VIP 2 Solves this by bringing the users hip joint closer to the reclining pivot point of the chair, which is located close to the human hip joint. It synchronises the motion of the body and the chair, allowing for minimal sliding.

Foldable and Durable

The VIP2 is light, foldable and easy to fit into a car boot. The additional push-bar reinforces the folding structure of the wheelchair.

Tension Adjustable Backrest

Unlike conventional wheelchairs, the backrest uses adjustable velcro straps to change the tension and conform according to the wheelchair user. This helps prevent slouching and the hammock effect.

J Hook Handles

Providing extra height and allowing the caregivers to move the VIP2 with ease and care.

Other Features

- Tension adjustable backrest
- Curved headrest included as standard
- Ergonomic, adjustable and detachable armrests
- Detachable swing-away footrest
- After-market backrest systems can be attached
- Microshield inhibiting bacterial and fungal growth
- Anti-Tippers
- Transit option may be available by special request



The center of gravity shifts between front and rear wheels.
• indicates COG of the chair

• indicates human hip joint
• indicates the reclining pivot point of the chair

+ ACCESSORIES

- Biomechanical elevating legrest
- Cane holder
- I.V Pole Holder
- Hemi Armpad
- Oxygen Tank Holder
- Oxygen tank holder

+ SPECIFICATIONS



Seat widths
16" (41cm) or
18" (46 cm)




Seat depth
16" (41cm) or
18.5" (46cm)




Wheels
7" (18cm) - front
20" (51cm) - rear



Seat to floor height
19.5" or 18.5"

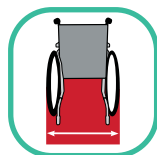

Overall width
26" (66cm) to
28" (71 cm)

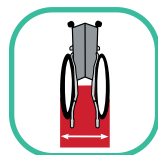

Folded width
38cm



Back height
49 / 54cm



Overall length
107-112cm



Overall height
138cm



Total weight
25.8kg



User weight limit
115kg