



meywalk® 4

The Dynamic Walker

18 months old to Adult



Quick & easy

3 Clicks and you're in!

Paediatric

Use at School, at Home, Indoors or Outdoors.

Adult

Adult Ideal in the Community, Acute settings, Residential/Care Homes for Falls Prevention, Single Handed Care, Stroke Rehab or Supporting Independence.



The Meywalk4 mimics natural gait and with its "3 clicks and you're in" easy transfer system, ensuring you get more time walking

- **Dynamic Pelvic Support:** Allows for natural hip rotation and 360° controlled steering
- **Sprung Loaded 'Springy' Frame:** Mimics natural gait
- **Easy Transfers:** 3 Clicks and you're away
- **Multi-Purpose:** Use it facing Backwards "Disco Position" for Close Quarter Therapy and Activity
- **Tilt In Space:** Posterior or Anterior Pelvic support
- **Tool Free Adjustment:** No more allen keys or spanners!



*2 clicks for Micro and Mini

Meywalk 4 Micro

Meywalk 4 can be used
'reversed' for close-up work



meywalk®4

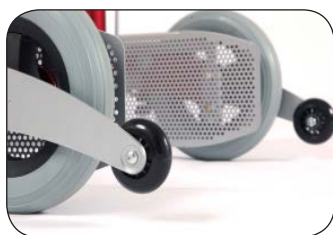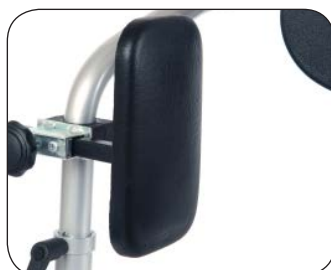








The right match

An extended range of accessories are available for Meywalk®4 ensuring the walker will match the specific needs of the user:

- **Trunk pad**
- **Hip pads with trunk pad**
- **Leg divider**
- **Leg guides**
- **Anti-tip support**

Technical data	Micro	Mini	Small	Medium	Large
Length:	77 cm	80 cm	85 cm	93 cm	110 cm
Outer width:	51 cm	56 cm	62 cm	62 cm	68 cm
Access width (inner width):	42 cm	44 cm	49 cm	50 cm	57 cm
*Seat height:	25-45 cm	28-52 cm	50-67 cm	60-84 cm	80-104 cm
*Trunk support height:	52-65 cm	55-77 cm	74-99 cm	101-133 cm	120-152 cm
Trunk support circumference	70, 90cm	70, 90, 105 cm	70, 90, 105, 120 cm	70, 90, 105, 120, 130 cm	70, 90, 105, 120, 130 cm
Max load:	60 kg	80 kg	100 kg	130 kg	130 kg
Weight:	13 kg	20 kg	25 kg	32 kg	34 kg
Colour:	Black	Yellow	Green	Red	Black
Order number:	WA1202	WA1200	WA1201	WA1203	WA1204

Seven spring user weights:

- | | | | |
|---|---|---|---|
|  Orange 10 - 15 kg |  Purple 15 - 20 kg |  White 20 - 50 kg |  Yellow 45 - 65 kg |
|  Blue 65 - 80 kg |  Red 80 - 100 kg |  Silver 100 - 130 kg | |