

# Pediatric Seating and Positioning System



As we recognize the important role of early intervention in cognitive development, more and more children use mobility bases at an earlier age. *VARILITE* JUNIOR seating and positioning products were designed with the knowledge that kids are not small adults—that's why we focused on the not-so-small details that make our products just for them.

The VARILITE JUNIOR line features a full-range of size-appropriate seating and positioning items to support children in their mobility device, allowing them to go where their imagination takes them. By exploring their world, children can gain confidence, increased peer interaction and independence. And VARILITE JUNIOR seating systems are there to support them along the way.



### **Sized For Children**

VARILITE JUNIOR products were designed with children in mind using anthropometric data.



## **Designed for Growth**

VARILITE JUNIOR products allow for growth with simple-to-adjust hardware.



## Lightweight

VARILITE JUNIOR products combine lightweight positioning materials to support kids on the move.



### Low Maintenance for Busy Parents

VARILITE JUNIOR products are simple to use to keep up with busy families.

## **VARILITE JUNIOR Seat Cushion**

The JUNIOR Seat Cushion provides a comfortable, stable foundation for exploration and adventure. Moderate skin protection and positioning is achieved using *VARILITE* air-foam floatation technology and a soft tailored cover.

Part #	Description	Dimensions
77000110102	JUNIOR 10x10 seat cushion w/ incontinence cover	10x10 in. (25x25 cm)
77000110122	JUNIOR 10x12 seat cushion w/ incontinence cover	10x12 in. (25x30 cm)
77000111112	JUNIOR 11x11 seat cushion w/ incontinence cover	11x11 in. (28x28 cm)
77000111132	JUNIOR 11x13 seat cushion w/ incontinence cover	
77000112122	JUNIOR 12x12 seat cushion w/ incontinence cover	12x12 in. (30x30 cm)

**NOTE:** Height is 2.9 in. (7.4 cm) and weight is approximately 9.2 oz (261g) for all cushion sizes.

### Features



**Multi-stiffness foam** and a non-intrusive medial thigh separator promotes immersion, stability and positioning.



Water resistant, washable covers and a waterproof inner cushion let kids be kids.



**Quick adjustments** can be made using the self-inflating valve.

# **VARILITE JUNIOR Back Support**

Lightweight and easy to use, the JUNIOR Back Support is perfect for children and families on the go. Moderate positioning is achieved using an adjustable *VARILITE* air-foam floatation cushion and versatile mounting hardware. Offered in mid and deep contoured styles.

#### **Features**



Compact, minimal design directs focus to the child, not the equipment.




**Simple adjustments** allow the back system to grow with your child.





**Lightweight** air-foam construction coupled with lightweight aluminum allows your child to manually propel with ease.



**Simple operation** means less time fussing with hardware and more time on the go.

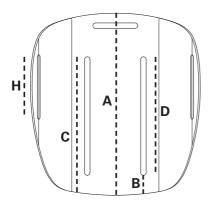

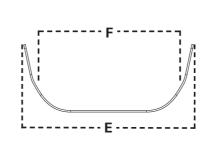

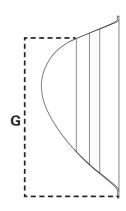

	Description	Mid Back	Deep Back
A	Back height*	9.25 in. (23.5 cm)	9.25 in. (23.5 cm)
В	Bottom of back to beginning of vertical slot	1.0 in. (2.5 cm)	1.0 in. (2.5 cm)
С	Bottom of back to top of vertical slot	7.0 in. (17.8 cm)	7.0 in. (17.8 cm)
D	Range of hardware adjustment in vertical slot	4.91 in. (12.5 cm)	4.91 in. (12.5 cm)
Е	Outside width of shell	8.6 in. (21.8 cm)	9.9 in. (25.1 cm)
F	Inside width of support surface	7.6 in. (19 cm)	8.9 in. (22.6 cm)
G	Height of fixed lateral support	N/A	5.7 in. (14.5 cm)
Н	Height of lateral slot	3 in. (7.6 cm)	N/A
Fits wheelchairs sizes		9-11 in. (23-28 cm)	9-11 in. (23-28 cm)
Weight, including hardware		2.4 lbs. (1.1 kg)	2.9 lbs. (1.3 kg)
Part #		4430111020	4430411020

<sup>\*</sup>Shell only. Complete back system height may vary by 1 in. (2.5 cm) due to the cushion rollover.

©2019 Cascade Designs, Inc.

4000 1st Ave. S, Seattle, WA 98134 Phone: (800) 827-4548 Fax: (206) 343-5795 **Authorized EU Representative:** 

Cascade Designs, Limited
Dwyer Road, Midleton, Co. Cork
Republic of Ireland



Document # 02331-2

