

Move More. Live More.

Experience the world with WHILL and live life without boundaries.

Innovative products built to challenge the perception of what a mobility device can be.

WHILL Model C2





Convenience

The Model C2 is simple to operate and can make tight turns to get you wherever you need to go.



The basket is conveniently located under the seat, making it easy to get your items in and out.

Dimensions: 24 cm x 29 cm x 43 cm



Lift the arm to easily sit down in the chair from the side. Stand on the foot support for added stability.



Adventure with WHILL - the Model C2 can make tight turns, allowing you to go places that a typical scooter can't take you.



Turn the power on and adjust the speed with one hand. The device will automatically brake when you release the controller, even on an incline.






Powerful

Impressive power that allows you to see and do more. Live the life you aspire to with the WHILL Model C2.

Get out and adventure!



Effortlessly climb over obstacles up to 5 cm in height with two powerful motors and large front omni-wheels.

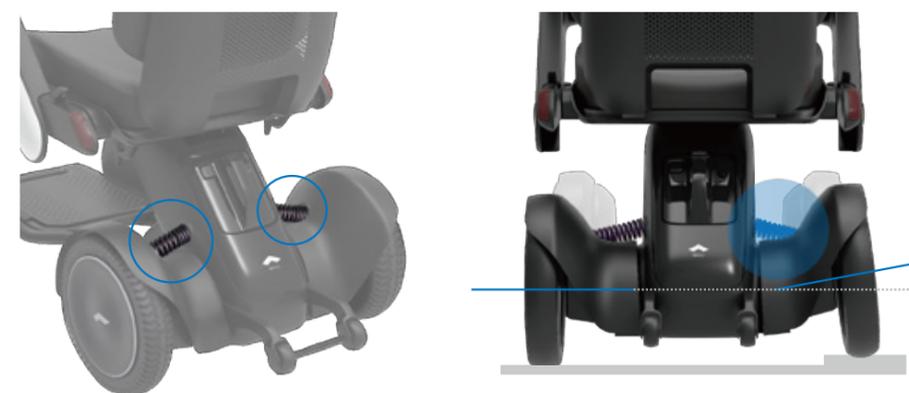

Rear wheel suspension absorbs the impact when driving over bumpy roads or obstacles.



Easily disassemble
to fit in the trunk
of your car

*Image shown above are for
illustrative purposes only


Portable

Disassembles in four easy steps.

The Model C2 can be disassembled for easy transportation in the trunk of any car. Whether you're shopping on the weekend or traveling with your family, the Model C2 allows you to live life to the fullest.



1 Remove the battery.

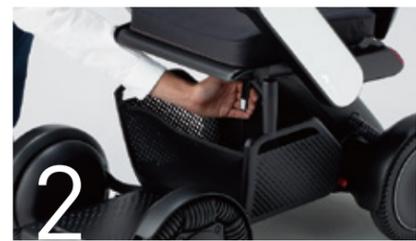

2 Raise the lever under the seat.



3 Pull the seat upwards.



4 Pull the lever on the side of the battery toward you to separate the parts.



We've thought of it all!

Feedback from our customers has helped shape the Model C2. We've packed in all the features you're looking for, including a remote lock and the ability to adjust the max speed for added safety.



Tail Light

The tail light is clearly visible even with a bag hung over the back support, keeping you safe on the road at night.

Prevents drifting



Drive straight on a tilted path.

The controller can be positioned on the left or right side.



Battery
Detachable battery.

5

hour charge



Model C2 can travel 18 km on a 5-hour charge, giving you ample power to go about your day.




2.7 kg

The battery can be removed and charged in a regular household outlet.




App

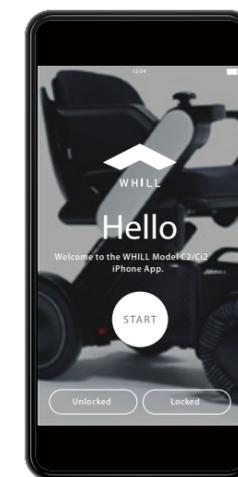


Check your total mileage, distance traveled, battery level, full charge cycles and more with the WHILL app*. You can also lock and unlock the unit.

*You will need to install the application beforehand.

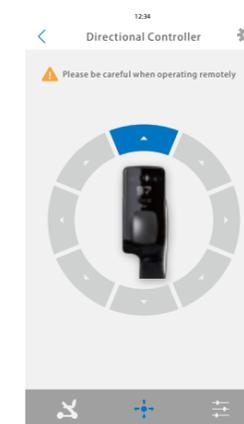


The intuitive controller screen allows you to control the unit remotely.

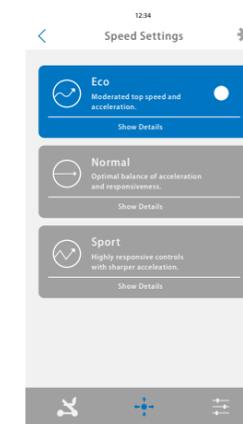

Set the drive mode to "Eco", "Normal" or "Sport" depending on your needs.

The Right Fit

The Model C2 can be easily adjusted from seat height and back support angle to controller position (right or left drive) to create a comfortable driving experience.

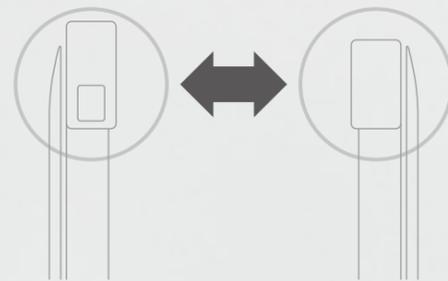

The controller can be positioned on the left or right side.

- 1. Controller Height 670, 690, 710, 730 mm
- 2. Back Support Angle 90, 95, 100°
- 3. Back Support Height ... 340 - 460 mm
- 4. Seat Height 405, 425, 445, 465, 485 mm
- 5. Seat Angle 5°
- 6. Seat Depth 400 mm
- 7. Seat Width 400 mm
- 8. Foot Support Angle 2°
- 9. Foot Support Height 127 mm

Warning: When using your own cushion, for safety, be sure to use one that has sufficient flame retardancy.

Performance Specifications

Driving Range 18 km
 Maximum Speed 6 km/h
 Drive System Two Motor Drive System
 Incline Capability 10°
 Turning Radius 760 mm
 Obstacle Clearance 50 mm
 Ground Clearance 57 mm
 Braking System Electromagnetic
 Max Weight Capacity 136 kg

Other Specifications
 Charging Time 5 hours
 Battery Lithium-ion (25.3 V x 10.6 Ah)
 Battery Charge Cycles ~1.000
 Battery Wh 269 Wh
 Operating Conditions -15 - 40°C
 Seat Cushions 400 mm x 400 mm
 450 mm x 450 mm
 Tail Lamps LED light (red)
 Horn Built-in, sound level > 65dB
 Wheelchair class A

Total length (including foot rest): 985 mm
 Total width: 554 mm
 Total height: 745 - 945 mm

Seat assembly

Length: 570 mm
 Width: 550, 600, 650 mm
 Height: 480 mm

Rear drive base

Length: 590 mm
 Width: 554 mm
 Height: 300 mm

Front drive base

Length: 720 mm
 Width: 554 mm
 Height: 270 mm

Type of front wheels: Omni-Wheel

Type of rear wheels: Standard: Airless Tire



*Do not use the device as a seat in an automobile or other vehicle.

*The unoccupied device is suitable for land and air transport.

Accessories

Purchase accessories for even more convenience and personalization.



Arm Covers

You can change the color by purchasing an additional one from the eight colors.

Left and right pair.



Lap Belt

Firmly stabilizes the sitting position even on slopes.



Battery

The 5.9 lbs lightweight battery can also be used as a spare battery.



Battery Charger

A great compliment to a spare battery, it helps let you go from one fully charged battery to the next fully charged battery to keep you moving.



Dust Cover

Protects the device from dust and dirt when stored outdoors*.

*It is not waterproof. Do not use for protection from rain.



Cane Holder

T-shaped canes, clutch-type canes, and crutches can be fixed to the main body.



Joystick

You can operate it by grasping it with your hand or sandwiching it between your fingers.



Side Bag

A side bag that can be attached to the arm support. You can store important things such as smartphones and wallets near your body. There is also a zipper pocket inside the bag.



Armrest Pad

A pad with cushioning performance that can be attached to the arm support. The height of the arm support can be raised by 3 cm.



Move More.
Live More.

